

WelchAllyn®



Welch Allyn  
CP 150™ Electrocardiograph  
Simple. Fast. Connected.



# Help your department run more efficiently with the Welch Allyn CP 150 Resting ECG—now with optional wireless capability.

Simple, streamlined ECG testing designed to help improve clinician workflows and patient care, providing comprehensive ECG data at the right time.

The CP 150 combines the best traditional box ECG features (e.g., full-size printouts and clinical decision support with MEANS interpretation) with connectivity enhancements designed especially for today's fast-paced environments.

In fact, **94% of CP 150 users surveyed** say the touchscreen display and simple user interface have helped them improve their efficiency.\*

\* Source: Survey by Walker Information, Inc. of CP 150 users.

**94%**  
CP 150 users surveyed said:  
"I work more efficiently with CP 150"

## Simple



### Streamline workflow

Can help improve accuracy, reduce errors and streamline workflow with automatic patient identification and worklists.



### Intuitive touchscreen display

User-friendly display and interface helps minimise training and save time.



### Lead quality display

Help avoid repeat ECGs by identifying lead connection problems before completing a test.

## Fast



### Instant-On Feature

Help eliminate typical start-up delays.




### DICOM® or HL7® Worklists

Can help reduce patient data re-entry by importing the DICOM or HL7 worklists.



### 12-Lead ECG Preview

Large display provides three, six or 12-lead ECG preview options.

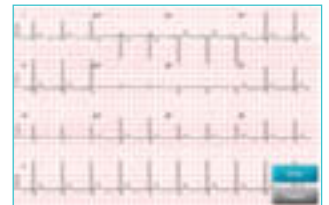

### ECG Print Preview

Can help save time and reduce paper costs by verifying the ECG test quality prior to printout.



LAN or Wi-Fi®

**Send ECG results to:**

-  NETWORK  
Network folder as PDF  
(For import into EMR)
-  HL7  
HL7 server  
(EMR, CVIS)
-  DICOM  
DICOM image server  
(EMR, CVIS, PACS)

## Connected

### LAN or optional wireless connectivity aimed at increasing device mobility and utilisation

- Optional wireless connectivity through an integrated 802.11 radio
- Share results immediately to put the data where you need it—in your EMR/CVIS/PACS—right away
- Meets globally accepted data encryption standards for enhanced patient information security

### Connectivity options tailored to enhance your practice's workflow efficiency

#### 1. Wireless PDF transfer to network folders

#### 2. Compatible with HL7 standards

- Bidirectional workflow with your EMR/CVIS<sup>1</sup>
- Help reduce patient data re-entry by importing the HL7 worklist
- Accepts standard ADT (Admit\Transfer\Discharge) messages as well as ORU (Observational Report Unsolicited) and ORM (General Order Message)

#### 3. DICOM® conformance can enhance interoperability with multiple DICOM-compatible platforms

- Bidirectional workflow with your EMR/CVIS<sup>1</sup> or PACS<sup>2</sup> without the middleware
- Help reduce patient data re-entry by importing the DICOM modality worklist
- Upgrade your existing CP 150 to make it DICOM compatible

*HL7 is the registered trademark of Health Level Seven International. DICOM is the registered trademark of the National Electrical Manufacturers Association for its standards publications relating to digital communications of medical information. Wi-Fi is a registered trademark of Wi-Fi Alliance.*

1. Cardiovascular information system, 2. Picture archiving and communication system



CP 150 mounted on hospital cart with cable arm and shelf



Office cart wheels

## The optional mobile stand makes the CP 150 ECG easy to move from room to room.

- Two choices of mobile carts: Hospital Cart—with four durable rubber wheels, two fixed with brake or four free-wheeling, or Office Cart—with four free-wheeling plastic wheels
- Two side baskets for additional storage
- Optional cable arm and shelf to keep patient cable and printed tests organised

## Ordering Information

### Selected Device Configurations\*:

CP150-2EN4	CP 150 Non-Interpretive Electrocardiograph
CP150A-2EN4	CP 150 Interpretive Electrocardiograph
CP150W-2EN4	CP 150 Non-Interpretive Electrocardiograph with Wireless Capability
CP150AW-2EN4	CP 150 Interpretive Electrocardiograph with DICOM
CP150D-2EN4	CP 150 Non-Interpretive Electrocardiograph with DICOM
CP150WD-2EN4	CP 150 Non-Interpretive Electrocardiograph with Wireless Capability and DICOM
CP150AWD-2EN4	CP 150 Interpretive Electrocardiograph with Wireless Capability and DICOM

### Accessories:

719654	10-lead Patient Cable, IEC, Banana, CP 50 / CP 150
721329	10-lead Patient Cable, IEC, Banana, CP 50 / CP 150 - 1,5 Meter
721361	10-lead patient cable, IEC, Pinch, CP 50 / CP 150 (1 m trunk length)
45008-0000	Resting Tab Electrodes, 1000/Box
715006	ECG Multifunction Electrode Adaptors, 10/Pack
105353	Recording Paper, CP 100/150/200, 200 Sheets/Pack, 5 Packs/Case
BATT99	9-Cell Lithium Ion Battery, Single Pack

### Accessories—Mobile Carts:

105341	CP 150 Office Cart
105342	CP 150 Hospital Cart
105343	CP 150 Cable Arm & Shelf Cart

### Welch Allyn Partners In Care<sup>SM</sup> Extended Warranty Programs:

S4-CP150	CP150, Warranty Extension 1 year
S4-CP150-2	CP150, Warranty Extension 2 years
S4-CP150-5	CP150, Warranty Extension 5 years

### Upgrades:

105410	Interpretation Upgrade, CP 150
406814	CP 50/150 Connectivity Kit
106736	CP 150 DICOM Upgrade Kit

### Warranty Information:

CP 150 Electrocardiograph—3 years  
Patient Cable/Battery—90 days

Contact your Welch Allyn representative or visit [www.welchallyn.co.uk](http://www.welchallyn.co.uk) for more information.

Welch Allyn UK Ltd.  
Cublington Road  
Aston Abbotts  
Buckinghamshire  
HP22 4ND  
(Tel:) 0207 365 6780  
(Fax) 0207 365 9694

[WWW.WELCHALLYN.CO.UK](http://WWW.WELCHALLYN.CO.UK)

© 2017 Welch Allyn SM3052UKI RevC MC13651 2017-01

**WelchAllyn**<sup>®</sup>

**Hill-Rom**<sup>®</sup>